

# BEAKER DYER UP TO 135°C AND WASHTESTER

## TYPE ALT-II

### GENERAL DESCRIPTION

Laboratory apparatus for high temperature exhaust dyeing up to 135°C. Textile samples are kept in stainless steel beakers made of one piece only (no welding). Heating by means of electric infrared heaters.


Automatic dye process control, with Pico Touch screen controller to programm heating and cooling gradients, holding times, alarms (memory for 99 programs with 20 steps each). Optional: memory for 30 programs.

Optional: Dosing during dye process.

Wastests according to norms beside dyeing possible.

### TECHNICAL DATA

Heating	5,0 kW
Rotation of beakers	30 - 60 rpm
Temperature gradient	0,3 - 4,0 C/min
Dimensions	L 83 x P 70 x A 94 cm
Weight	165 kg

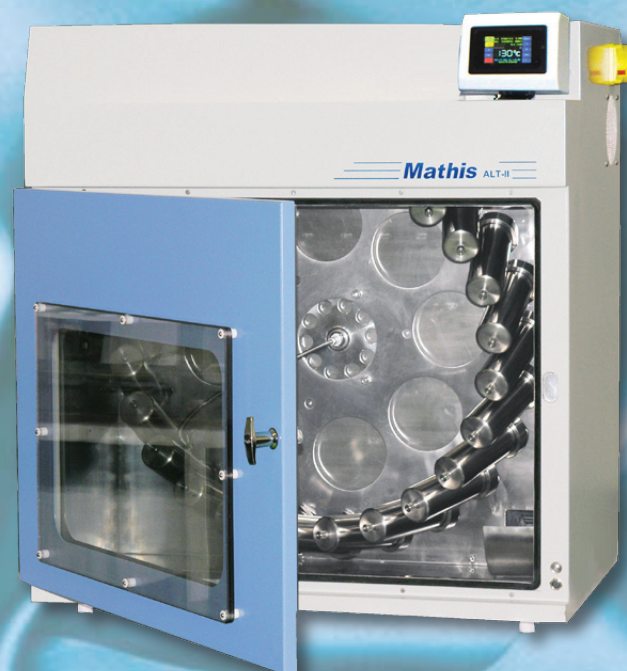

Optional:  
Dosing during process



Dosing syringe  
5,0 ml

## ACCESSORIES ALT-II:



holder with 32 x 100ml beakers  
with or without dosing nipple




holder with 24 x 200/300ml beakers  
with or without dosing nipple




holder with 8 x 500ml beakers  
for wash tests



holder with 12 x 500ml beakers  
with or without dosing nipple

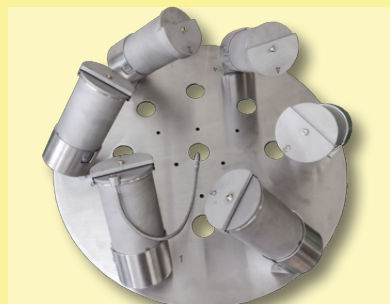

holder with 6 x 1000/1200ml  
or 1400ml beakers  
with or without dosing nipple

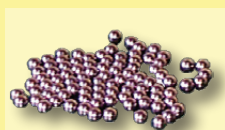

stainless steel balls  
for wash tests


yarn holders  
for 100, 200, 300, 500,  
1000, 1200, 1400 or  
1800 ml beakers



holder with 6 x 1800ml beakers  
with or without dosing nipple



holder with 1 x 5 liter beaker  
with or without dosing nipple

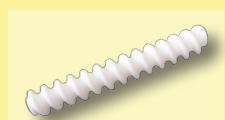

special liquor  
displacement body  
for short liquor  
rations

### Werner Mathis AG

Rütisbergstrasse 3  
CH-8156 Oberhasli/Zürich, Switzerland  
Telefon 41(0)44-852 5050  
Telefax 41(0)44-850 6707  
E-Mail info@mathisag.com  
Homepage www.mathisag.com

### Aparelhos de Laboratório Mathis Ltda.

Estrada dos Estudantes, 261  
06707-050 Cotia - SP, Brasil  
Tel: 55 (0)11 3883-6777  
Fax: 55 (0)11 3883-6750  
E-mail: mathis@mathis.com.br  
Site: www.mathis.com.br